



Performance Check of Turbidity Meter

Equipment Ref. No. : ET/0505/010 Manufacturer : HACH
Model No. : 2100Q Serial No. : 11110 C 014260
Date of Calibration : 07/01/2014 Due Date : 06/04/2014

Theoretical Value of Turbidity Standard (NTU)	Measured Value (NTU)	Difference % *
20	19.2	-4.08
100	104	3.92
800	793	-0.88

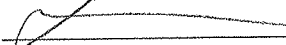
(*) Difference = (Measured Value – Theoretical Value) / Theoretical Value

Acceptance Criteria

Difference : -5 % to 5 %

The turbidity meter complies * / ~~does not comply~~ * with the specified requirements and is deemed acceptable * / ~~unacceptable~~ * for use. Measurements are traceable to national standards.

Prepared by : 

Checked by : 



Internal Calibration Report of Dissolved Oxygen Meter

Equipment Ref. No. : <u>ET/EW/008/004</u>	Manufacturer : <u>YSI</u>
Model No. : <u>Pro 2030</u>	Serial No. : <u>10F 101978</u>
Date of Calibration : <u>29/01/2014</u>	Calibration Due Date : <u>28/04/2014</u>

Temperature Verification

Ref. No. of Reference Thermometer : ET/0521/008

Ref. No. of Water Bath : ---

		Temperature (°C)		
Reference Thermometer reading	Measured	20.2	Corrected	19.8
DO Meter reading	Measured	19.6	Difference	0.2

Standardization of sodium thiosulphate (Na₂S₂O₃) solution

Reagent No. of Na ₂ S ₂ O ₃ titrant	CPE/012/4.5/001/8	Reagent No. of 0.025N K ₂ Cr ₂ O ₇	CPE/012/4.4/001/24
		Trial 1	Trial 2
Initial Vol. of Na ₂ S ₂ O ₃ (ml)		0.00	10.50
Final Vol. of Na ₂ S ₂ O ₃ (ml)		10.50	20.95
Vol. of Na ₂ S ₂ O ₃ used (ml)		10.50	10.45
Normality of Na ₂ S ₂ O ₃ solution (N)		0.02381	0.02392
Average Normality (N) of Na ₂ S ₂ O ₃ solution (N)		0.02387	
Acceptance criteria, Deviation		Less than ± 0.001N	
Calculation: Normality of Na ₂ S ₂ O ₃ , N = 0.25 / ml Na ₂ S ₂ O ₃ used			

Linearity Checking

Determination of dissolved oxygen content by Winkler Titration *

Purging Time (min)	2		5		10	
	1	2	1	2	1	2
Trial						
Initial Vol. of Na ₂ S ₂ O ₃ (ml)	0.00	11.30	22.40	0.00	8.20	13.20
Final Vol. of Na ₂ S ₂ O ₃ (ml)	11.30	22.40	30.80	8.20	13.20	18.10
Vol. (V) of Na ₂ S ₂ O ₃ used (ml)	11.30	11.10	8.40	8.20	5.00	4.90
Dissolved Oxygen (DO), mg/L	7.24	7.11	5.38	5.25	3.20	3.14
Acceptance criteria, Deviation	Less than + 0.3mg/L		Less than + 0.3mg/L		Less than + 0.3mg/L	

Calculation: DO (mg/L) = V x N x 8000/298

Purging time, min	DO meter reading, mg/L			Winkler Titration result *, mg/L			Difference (%) of DO Content
	1	2	Average	1	2	Average	
2	7.23	7.11	7.17	7.24	7.11	7.18	0.14
5	5.33	5.17	5.25	5.38	5.25	5.32	1.32
10	3.30	3.12	3.21	3.20	3.14	3.17	1.25
Linear regression coefficient						0.9990	



Internal Calibration Report of Dissolved Oxygen Meter

Zero Point Checking

DO meter reading, mg/L	0.00
------------------------	------

Salinity Checking

Reagent No. of NaCl (10ppt)	CPE/012/4.7/002/15	Reagent No. of NaCl (30ppt)	CPE/012/4.8/002/15
-----------------------------	--------------------	-----------------------------	--------------------

*Determination of dissolved oxygen content by Winkler Titration ***

Salinity (ppt)	10		30	
Trial	1	2	1	2
Initial Vol. of Na ₂ S ₂ O ₃ (ml)	0.00	12.00	24.30	35.30
Final Vol. of Na ₂ S ₂ O ₃ (ml)	12.00	24.30	35.30	46.50
Vol. (V) of Na ₂ S ₂ O ₃ used (ml)	12.00	12.30	11.00	11.20
Dissolved Oxygen (DO), mg/L	7.69	7.88	7.05	7.18
Acceptance criteria, Deviation	Less than + 0.3mg/L		Less than + 0.3mg/L	

Calculation: $DO (mg/L) = V \times N \times 8000/298$

Salinity (ppt)	DO meter reading, mg/L			Winkler Titration result**, mg/L			Difference (%) of DO Content
	1	2	Average	1	2	Average	
10	7.81	7.92	7.87	7.69	7.88	7.79	1.02
30	7.33	7.42	7.38	7.05	7.18	7.12	3.59

Acceptance Criteria

- (1) Differenc between temperature readings from temperature sensor of DO probe and reference thermometer : < 0.5 °C
- (2) Linear regression coefficient : >0.99
- (3) Zero checking: 0.0mg/L
- (4) Difference (%) of DO content from the meter reading and by winkler titration : within ± 5%

The equipment complies # / does-not-comply # with the specified requirements and is deemed acceptable # / unaccepttable # for use.

Delete as appropriate

Calibrated by

:

Approved by:



Performance Check of Salinity Meter

Equipment Ref. No. : ET/EW/008/004 Manufacturer : YSI
Model No. : Pro 2030 Serial No. : 10F 101978
Date of Calibration : 29/01/2014 Due Date : 28/04/2014

Ref. No. of Salinity Standard used (30ppt)

S/001/5

Salinity Standard (ppt)	Measured Salinity (ppt)	Difference %
30.0	32.3	7.38

Acceptance Criteria

Difference : <10 %

The salinity meter complies * / ~~does not comply~~ * with the specified requirements and is deemed acceptable * / ~~unacceptable~~ * for use. Measurements are traceable to national standards.

Checked by : 

Approved by : 